

Instructions for use

Please read this carefully because it will help you to use HYLO DUAL INTENSE® correctly.

HYLO DUAL INTENSE®

0.2 % Sodium hyaluronate and 2 % ectoin

For intensive therapeutic lubrication of dry eyes and stabilisation of the tear film.

Relieves inflammatory symptoms such as irritation, itching and burning.

PRODUCT DESCRIPTION

HYLO DUAL INTENSE® are aqueous, sterile, preservative-free, phosphate-free eye drops.

COMPOSITION

1 ml of HYLO DUAL INTENSE® contains 0.2 % sodium hyaluronate, 2 % ectoin, boric acid, borax and water for injections.

INTENDED PURPOSE

For the intensive therapeutic lubrication of dry eyes and stabilisation of the tear film. Relieves inflammatory symptoms such as irritation, itching and burning.

MEDICAL INDICATION

Dry eye disease

SUITABLE FOR

Can be used by adults and children older than 2 years. Please see special dosage recommended in section "Dosage" below.

Another person may assist in the application of the product if the patient is a child, an elderly person or a person who is unable to use the product themselves (e.g., due to injury, illness or limited mobility).

Based on current knowledge, the product can be used during pregnancy or lactation. Please consult your eye doctor before use and see special dosage recommended in section "Dosage" below if you are pregnant or breastfeeding.

CLINICAL BENEFIT, PERFORMANCE AND EFFECTIVENESS

Applying HYLO DUAL INTENSE® provides an effective and long-lasting lubrication of the eye surface. This will relieve the inflammatory symptoms of dry eye disease, including prolonged or more severe eye irritation and dryness, foreign-body sensation, burning or tearing eyes.

HYLO DUAL INTENSE® contains sodium hyaluronate and ectoin. Sodium hyaluronate is a natural substance that forms a long-lasting lubricating film on the surface of the eye. Ectoin, a natural substance with water-binding properties, stabilizes the lipid layer of the tear film. This protects against excessive evaporation of tears. The high viscosity and dual action of HYLO DUAL INTENSE® provides long-lasting relief for severe and persistent symptoms.

CONTRAINDICATIONS

Do not use if you are hypersensitive to any of the ingredients.

Do not use in children younger than 2 years.

It is advised not to use while wearing contact lenses due to the high viscosity of the eye drops.

SIDE EFFECTS

The product is preservative-free and well tolerated even after prolonged use. In rare cases, adverse reactions like hypersensitivity, burning or excessive tearing, pain, foreign body sensation or other transient local irritations have been reported which resolved when the product was discontinued.

Temporary blurred vision may occur immediately after application due to the high viscosity of the eye drops. In rare cases this may cause a temporary headache in sensitive patients.

If any serious incident occurs during or as a result of the use of this product, please contact the manufacturer () and your competent authority.

PRECAUTIONS AND WARNINGS

If your symptoms persist or get worse after more than a week of using the product, consult your eye doctor or pharmacist.

Take care when using HYLO DUAL INTENSE® simultaneously with other eye drops.

Apply eye drops first and wait at least 30 minutes before using HYLO DUAL INTENSE®. HYLO DUAL INTENSE® should be used last. Eye ointments should, however, always be administered after using HYLO DUAL INTENSE®. Wait at least 30 minutes between applications.

HYLO DUAL INTENSE® contains boric acid and sodium borate in small amounts that are safe if the product is used correctly according to the instructions for use.

One bottle should only be used by a single person ().

DOSAGE

Use the product as needed, depending on the severity of your individual symptoms.

Generally, adults and children older than 2 years can apply one drop in each eye three times a day. If you need to use the product more frequently (e.g. more than 10 times per day), or if your symptoms persist for more than a week or get worse, see your eye doctor or pharmacist.

Because the product contains borate buffer, pregnant and breast-feeding women, and children older than 2 years should not use HYLO DUAL INTENSE® more than 10 times per day in each eye.

Suitable for long term use.


GENERAL INSTRUCTIONS FOR USE

Do not use the product if any part of the seal on the folding carton is broken or damaged before first use ().


Only use the product if it is intact.

The product is filled in a special COMOD bottle, a multi-dose container with an airless pump designed to keep the eye drops sterile without using preservatives.


The handling is shown in the following:

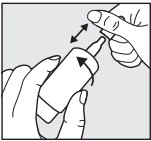


Step 1:
Remove the cap before each use.
For this purpose, hold the cap that protects the dropper tip right above the bottle. With your other hand, twist the bottle away from the cap.
Note: Before the first application, the pump mechanism must be activated. Turn the bottle upside down and press repeatedly on the bottom of the bottle until the first drop appears at the dropper tip. The bottle is now ready for further use.

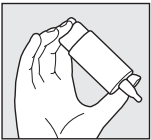


Step 2:
Hold the bottle upside down with the thumb on the shoulder of the bottle and the other fingers on the bottom of the bottle.

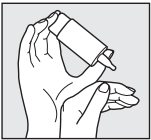

Step 3:
Support the hand holding the COMOD bottle with the other hand.

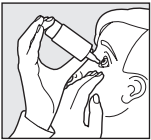

Step 4:
When applying the eye drops, avoid any contact between the dropper tip and your eye or skin.
Tilt your head back slightly and gently pull down your lower eyelid with your free hand. Press firmly and quickly on the bottom of the bottle. This activates the release mechanism to dispense a single drop. Due to the special COMOD valve system, the size and speed of each drop remain constant regardless of the pressure applied to the bottom of the bottle. Close your eyes after application so that the solution is evenly distributed over the eye surface.

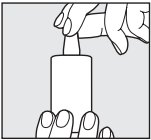

Step 5:
Always ensure that the dropper tip is dry and replace the cap immediately after use. This prevents residual liquid from crystallising and contaminating the dropper tip.

Contact the manufacturer (), your pharmacist or eye doctor if you have any questions about the correct use of this medical device.

10 ml of solution correspond to approximately 300 drops.
For technical reasons, a small amount of product remain in the bottle at the end of use.

STORAGE

Do not store above 25 °C ( 25 °C).
The product can be used for 6 months after first opening.
Do not use after the expiry date ().
Keep the product out of the sight and reach of children.
Dispose of this medical device with your household waste in accordance with local regulations. This will help to protect the environment.

THE FOLLOWING PACKAGE SIZES ARE AVAILABLE

- 10 ml of eye drops
- 7.5 ml o eye drops (sample pack)

DATE OF INFORMATION

July 2026
 URSAPHARM Arzneimittel GmbH
Industriestraße 35, 66129 Saarbrücken, Germany
www.ursapharm.de
 **And UK Responsible Person:**
Scope Ophthalmics Ltd.
1st Floor, Unit 4 Amberley Court
Whitworth Road, County Oak Way
Crawley, West Sussex, RH11 7XL
Freephone: 0808 175 3200
Fax: +44 203 515 9461
Email: info@scopeeyecare.com
www.scopeeyecare.com

0197





Symbol:	Explanation:
	Distributor
	Medical device
	CE mark
	Unique Device Identifier
	Sterilized using aseptic processing techniques
	Lot number
	Date of manufacture [YYYY-MM]
	Consult instructions for use
	Contains hazardous substances

